




Technical Data Sheet

Automatic filling and sealing machine MPP Uno

★ Product features

- Seals and fills 1 straw per cycle: ideal for handling small volumes
- Capacity: approx. 4000 straws/hour
- Exchangeable hopper, capacity: 2400 MiniStraws (0.25 ml) or 1200 medium straws (0.5 ml)
- Fast, tool-free conversion between 0.25 and 0.5 ml straws and vice versa (1-2 minutes)
- Ultrasonic sealing unit, ultrasonic generator in an external casing that can be installed at up to 2.5 m (ON/OFF switch directly on the machine)
- Vacuum unit with integrated reservoir and automatic readjustment function; vacuum level can be easily modified via buttons
- Can be operated in cold room or cabinet or at room temperature as stand-alone or combined with MiniJet printer (CombiSystem)
- Very compact device which takes up little space; machine dimensions: 460 x 350 x 410 mm (W x H x D) plus ultrasonic generator
- Optimal control by the user: filling process is performed right in front of the user
- Automated process: easy to use, no extensive operator training necessary
- Well-proven construction: high reliability and virtually no downtime
- Ideal conditions for hygiene and safety
- Power supply: 115 V/60 Hz or 230 V/50 Hz

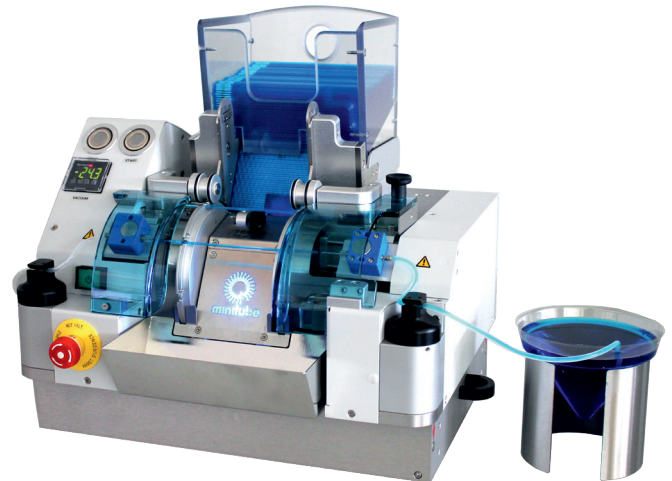

↗ Accessories

Straw hopper for 0.25 and 0.5 ml straws	13018/0006
Filling nozzle, for 0.25 ml straws	13017/0012
Suction nozzle, for 0.25 ml straws	13017/0022
Filling nozzle, for 0.5 ml straws	13017/0017
Suction nozzle, for 0.5 ml straws	13017/0027
Washer for 0.25 and 0.5 ml nozzles, 5000/bag	13021/3001
Filling and suction tubing, 380 mm, silicone, bevelled edge, 100/bag	13021/3800
Disposable semen cone, 100/bag	13018/0100

MPP Uno 13017/0000