




ePacker

Fully automated boar semen packaging machine

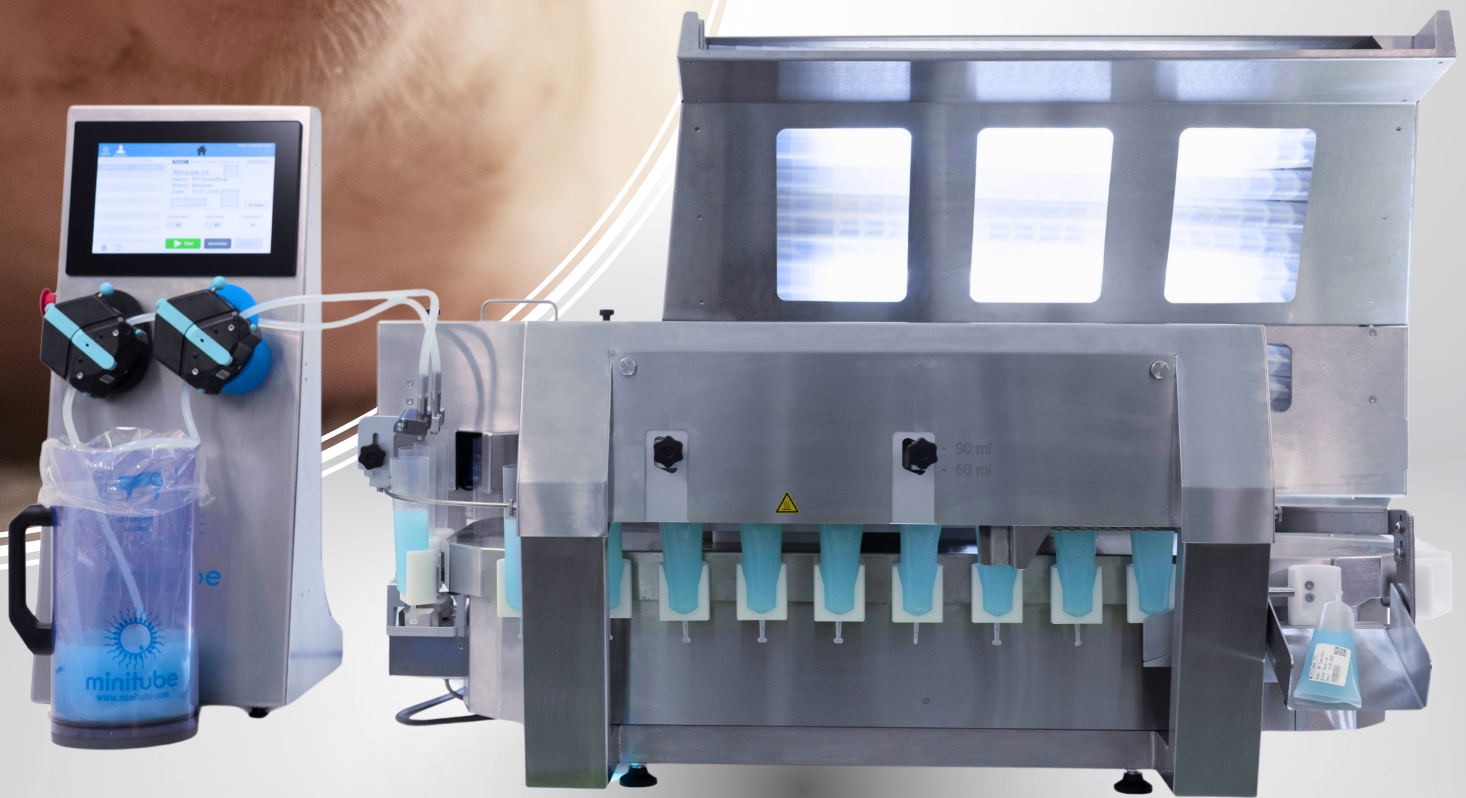

Robust, hygienic, quiet and fast, the ePacker is the perfect solution for boar studs looking for the latest technology!



High precision and speed

Equipped with two peristaltic pumps, the ePacker delivers a precise filling volume (+/- 1 ml) to each tube at a high processing speed of up to 1000 tubes per hour.



Autocalibration of fill volume

For expedited change between batches, minimizing downtime and maximizing effective production time.



Advanced sealing and printing

Thermal fusion technique (TFT) ensures secure and tight sealing. The integrated printing station provides state-of-the-art labels with text, logos, and codes.



Seamless and versatile operation

Versatile production with 60 ml and 90 ml tubes. The ePacker allows quick, tool-less adjustments for different tube sizes.



Quiet, user-friendly and hygienic design


Operating without compressed air, the ePacker is quieter than any other machine on the market. Its advanced automation, intuitive software and hygienic surfaces make it easy to operate, maintain and clean.



Economic value


The ePacker brings you the lowest cost of production, ensuring the highest profitability in your business.



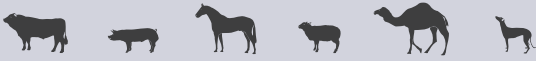

★ Product features

- Fills up to 1000 boar semen tubes (90 ml) per hour
- Equipped with two peristaltic pumps for a precise filling volume (+/- 1 ml)
- Auto-calibrates via load cells for consistent performance
- Uses a robust and novel thermal fusion technique (TFT) for secure and tight sealing
- Fills 60 and 90 ml tubes with short and flexible nozzle
- Allows quick, tool-less change between tube sizes in a few simple steps
- Modular design of key elements allows for easy access, maintenance and service
- Seamlessly connects to Prism10 and other lab software for data transfer of filling and labeling information, or operates stand-alone
- Supports a wide range of different label designs, including Barcode 128, QR, PDF417, DataMatrix, and more
- Holds up to 336 semen tubes (90 ml) in a high-capacity magazine with 16 stacks
- Features a hygienic design with stainless steel surfaces, rounded edges, and concealed cables for easy cleaning

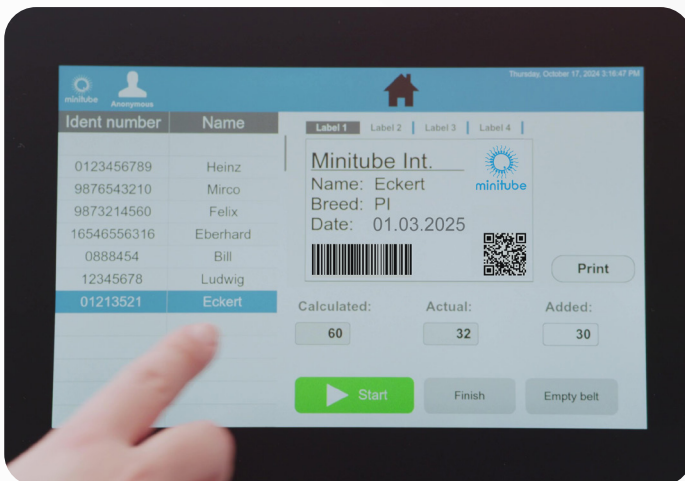

Unique WYSIWYG experience

What You See Is What You Get – The intuitive, touch-operated display provides a true real-time preview of the tube label, displaying the respective labeling information within the predefined label designs.

Order information

ePacker 13207/0000

Accessories

Dispensing tubing 13207/0160

Tube label, 50 x 32 mm, white, 2000/roll 13210/4964

